
## Canon

## Daily QA TOS-SS Phantom

Daily QA TOS-SS Phantom (180 mm)						
120 kV, 300 mA, 1 Sec.						
5.0 x 4, 10.0 mm (2 Stack)						
CFOV: Medium/ DFOV: 180 mm						
FC70, SU/FF/ <b>VFH</b>						
Measure using TOS-SS 1 Preset						
Record the Mean (HU) Value						

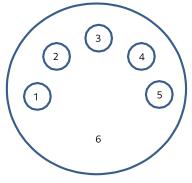

## AQ Prime SP TSX-303B

Normal Range (HU)					
1: Bone	910 to 970				
2: Fat	-115 to -95				
3: Air	-1020 to -980				
4: Tissue	115 to 135				
5: Oral	310 to 370				
6: Water	-5 to 5				

DATE	1 BONE	2 FAT	3 AIR	4 TISSUE	5 ORAL	6 WATER	TECH INITIALS
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							
27							
28							
29							
30							
31							